


Rapid CheckTM RSV Antigen Test

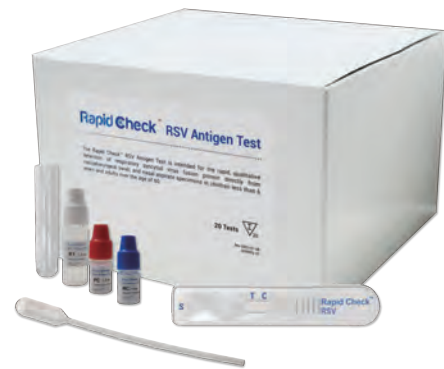

Product Description

A rapid, qualitative test for the detection of respiratory syncytial virus fusion protein directly from nasopharyngeal swab, and nasal aspirate specimens in children less than 6 years and adults over the age of 60.

The test is intended for use as an aid in the rapid laboratory diagnosis of acute respiratory syncytial virus infection in patients with symptoms consistent with RSV infection.

Swift Results, Seamless Care

Features/Benefits

Rapid Detection and Patient Management

RSV Infection results available in 15 minutes

Various Sample Types

Nasopharyngeal Swab and Nasal Aspirates from symptomatic patients less than 6 years and over the age of 60.

Simple and Effective

CLIA Waived testing at the point-of-care for use with minimally trained healthcare professionals

Premier Test Performance*

Positive Percent Agreement: 93.9%

Negative Percent Agreement: 97.7%

*Nasopharyngeal Swab and Nasal Aspirates. Refer to the Package Insert for additional performance claims.

Transport Medium Compatibility

Allows for analysis of preserved or transported specimens.

External Quality Controls Included

Facilitate routine quality inspection

Shelf Life that Spans Viral Seasonality


24 months from date of manufacture at room temperature storage

Rapid CheckTM RSV Antigen Test

Kit Contents

- 20 Individually packaged test devices
- 1 R1 (Extraction Reagent Solution)
- 20 Extraction Tubes
- 40 Disposable Pipettes
- 1 PC (Positive Control Reagent)
- 1 NC (Negative Control Reagent)
- 20 Sterile Nasopharyngeal Swabs
- 1 Instructions for Use

Test Procedure

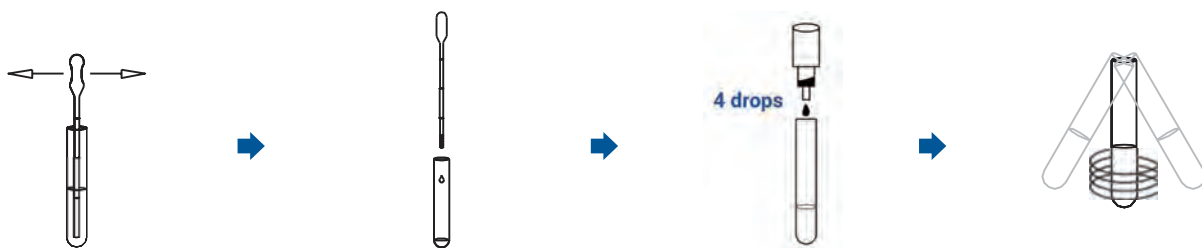


Draw 160 µL (second notch) of the sample up into a disposable pipette.

Add entire contents of the Pipette to the extraction tube.

Add four (4) **free falling** drops of R1 to the extraction tube while holding the bottle vertically over the tube.


Gently swirl the extraction tube to mix the sample and R1.

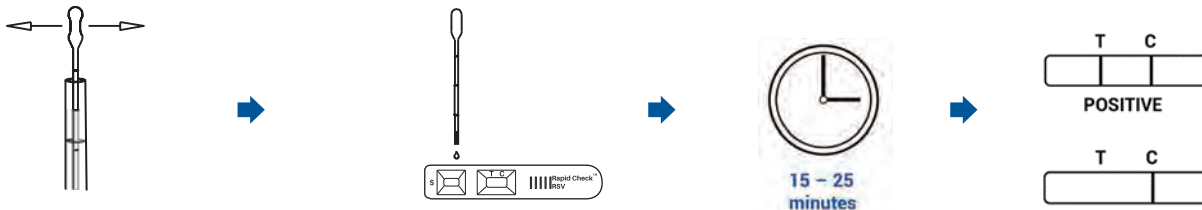

Draw 160 µL (second notch) of the sample up into a disposable pipette.

Squeeze entire contents of the pipette dropper onto the sample well of the test device.

Read test result between fifteen (15) and twenty-five (25) minutes after applying the extracted sample to the test device.

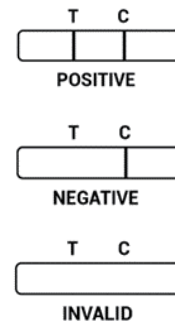

Ordering Information

REF	Description	QTY	UOM
GCRSV-502b-US	Test Kit	20T/Kit	Kit

***20 Sterile Nasopharyngeal Swabs included with kit**

www.germainelabs.com

GERMAINE[®]
LABORATORIES, INC.



800.854.8446
210.692.4192
210.692.4198 Fax